TAKE CARE!

SMART Baby Scale 2in1 with Bluetooth User Manual

GB	SMART Baby Scale 2in1 with Bluetooth	2
PL	Waga SMART 2w1 z Bluetooth	4
BG	Интелигентна електронна везна за бебета 2 в 1 с Bluetooth	6
CZ	Inteligentní elektronická kojenecká váha 2v1 s Bluetoothem	8
DE	Intelligente Babywaage 2-in-1 mit Bluetooth	10
EE	Nutikas elektrooniline imikukaal 2in1 Bluetooth-ga	12
ES	Báscula electrónica inteligente para bebés 2 en 1 con Bluetooth	14
FR	Le pèse-bébé électronique, balance connectée pour nourrisson 2en1 avec la technologie Bluetooth	16
GR	Έξυπνη ηλεκτρονική ζυγαριά για βρέφη 2 σε 1 με Bluetooth	18
HR	Inteligentna elektronička vaga za bebe 2 u 1 s Bluetoothom	20
HU	Okos, Buetoothszal rendelkező 2 az 1-ben babamérleg	22
IT	La bilancia elettronica per neonati 2 in1 dotata di Bluetooth	24
LT	Išmaniosios elektroninės svarstyklės kūdikiams 2w1 su Bluetooth įranga	26
LV	Inteligenti elektroniski zidaiņu svāri “divi vienā” ar Bluetooth tehnoloģiju	28
NL	Intelligente elektronische babyweegschaal 2in1 met Bluetooth	30
PT	Balança eletrónica inteligente 2em1 para bebés, equipada com Bluetooth	32
RO	Cântar electronic inteligent pentru nou-născuți 2in1 cu Bluetooth	34
RU	Умные электронные весы для новорожденных 2в1 с Bluetooth	36
SE	Smart elektronisk babyvåg 2i1 med Bluetooth	38
SK	Inteligentná elektronická váha na váženie dojčiat 2v1 s technológiou Bluetooth	40
UA	Розумні електронні ваги для новонароджених 2в1 з Bluetooth	42

Mum and baby logbook
Baby Fit

[Google Play](#)

[Apple Store](#)








To install the application, you may use a QR code or search for it using 'Baby Fit' name directly in the app store.

Before using this product, please read this manual.

INSTALLING THE APPLICATION

Scan the QR code or search for the app in the app store. Download it to your smartphone and open. Follow the instructions on the screen, enter your details and find the function of electronic scale. Perform the following steps:

1. Click to add a device.
2. Search and add the device, then click the displayed icon.
3. After successful pairing, go to the next step.
4. Enter the name of the device and select the location of its use - then start it.
5. Before the use, select the unit of measurement.
6. Select the user and enter the date.
7. Finish the measurement and save the data.

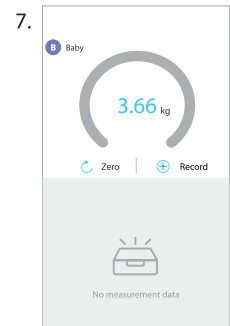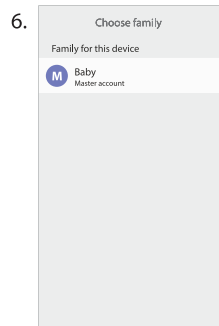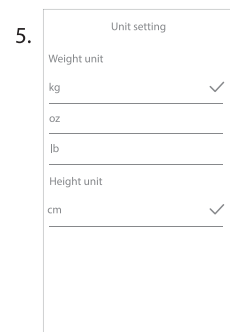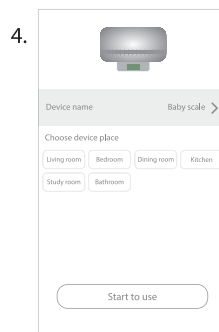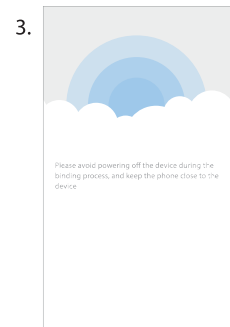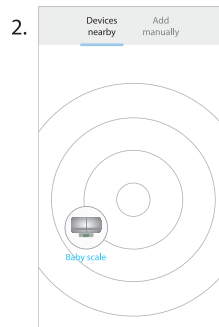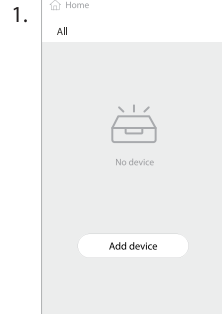PRODUCT FEATURES

- High accuracy sensor.
- White digits on a blue LCD display.
- Automatic shutdown.
- Overload warning.
- Function of units conversion.

TECHNICAL PARAMETERS AND MANUFACTURING STANDARDS

- Product Name: Electronic baby scale
- Model: NP - B1H
- Power supply: two AAA batteries
- Measuring range: 50 g - 100 kg
- Accuracy: 10 g
- Measurement error: + 3%
- Manufacturing standard: GB/T7722-2005


USAGE WAY

Batteries and their installation


The scale is designed as environmental-friendly and energy efficient device. It is supplied with two batteries. Before the use, remove foil from batteries and install them properly. If you see low battery message, replace batteries.

BATTERY REPLACEMENT:

- Open the battery compartment of the electronic scale.
- Insert new batteries place starting from one end of the battery and then push-in the other end (observe the polarity of batteries).

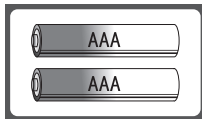DESCRIPTION OF BUTTONS

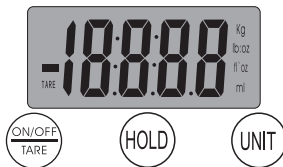A. Switch - the button is used for turning the device on and off. Press the switch to turn the device ON, and press it again to turn it OFF.

B. TARE - the button is used to tare weight after placing the tray on it (if supplied with the tray). Press 'Tare' button and then start weighing.

C. The device displays the weight in different units. To select a different unit, press UNIT.


The device allows you to record the results of measurements. After the weighing result is stable, press this button to save the result shown by the display. To delete the saved result, press the button again.

WEIGHING


1. Install the tray on the scale body


2. Turn on the scale by pressing ON / OFF button.

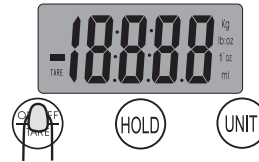3. When the display shows "0" value, place the baby on the scale and press "UNIT" button to change the unit of measure.

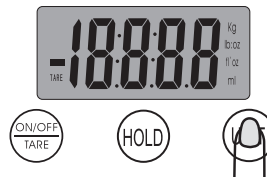4. Automatic shutdown of the scale (energy-saving mode).

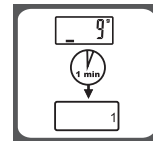If the scale does not detect any activity for one minute, it switches off automatically.

IMPORTANT INFORMATION

1. Before using the scale, place it on a flat and level surface
2. Do not use the scale in a place with a high humidity and very high or low temperatures.
3. Do not expose the scale to impacts and falls, as it contains electronic circuits of high precision.
4. For cleaning the scale, use only a soft cloth without chemicals. Do not flood the scale.
5. The scale is intended for domestic use only and can not be used for professional purposes.
6. If you do not plan to use the scale for a longer period, remove its batteries.
7. If the display shows "LO" message, immediately replace the batteries.
8. If the display shows "ERR" (overload), remove the weighted object from the scale, to prevent scale damage.

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